



ADVACAT

Purified Collagen (Chromic & Plain)
Natural Absorbable Surgical Suture

Consistent performance for
confidence in episiotomy repair























ADVACAT Purified Collagen (Chromic & Plain)

Natural Absorbable Surgical Suture

Consistent performance for confidence in episiotomy repair

Product Description

Code	Suture			Needle				Units Per Box
	Colour	Size	Length (cm)	Curvature	Type	Length (mm)	Speciality	
ADVA2003	YELLOW	2-0	152	-	-	-	-	12
ADVA2004	YELLOW	0	152	-	-	-	-	12
ADVA2005	YELLOW	1	152	-	-	-	-	12
ADVA2213	YELLOW	2-0	152	-	-	-	-	12
ADVA2214	BEIGE	0	152	-	-	-	-	12
ADVA2215	BEIGE	1	152	-	-	-	-	12
ADVA2216	BEIGE	2	152	-	-	-	-	12
ADVA4216	BEIGE	2-0	76	3/8 CIRCLE		30	-	12
ADVA4217	BEIGE	0	76	3/8 CIRCLE		30	-	12
ADVA4220	BEIGE	2-0	76	3/8 CIRCLE		40	-	12
ADVA4221	BEIGE	0	76	3/8 CIRCLE		40	-	12
ADVA4225	BEIGE	0	76	3/8 CIRCLE		45	-	12
ADVA4226	BEIGE	1	76	3/8 CIRCLE		45	-	12
ADVA4227	BEIGE	1	100	1/2 CIRCLE		45	HEAVY	12
ADVA4228	BEIGE	2	100	1/2 CIRCLE		45	HEAVY	12
ADVA4229	BEIGE	0	100	1/2 CIRCLE		45	HEAVY	12
ADVA4237	BEIGE	3-0	76	1/2 CIRCLE		20	-	12
ADVA4241	BEIGE	2-0	76	1/2 CIRCLE		30	-	12
ADVA4242	BEIGE	0	76	1/2 CIRCLE		30	-	12
ADVA4245	BEIGE	2-0	76	1/2 CIRCLE		40	-	12
ADVA4246	BEIGE	0	76	1/2 CIRCLE		40	-	12
ADVA4250	BEIGE	0	76	1/2 CIRCLE		45	-	12
ADVA4255	BEIGE	2-0	76	1/2 CIRCLE		30	HEAVY	12
ADVA4258	BEIGE	0	76	1/2 CIRCLE		40	HEAVY	12
ADVA4259	BEIGE	1	76	1/2 CIRCLE		40	HEAVY	12
ADVA4279	BEIGE	0	100	1/2 CIRCLE		40	HEAVY	12
ADVA5033	BEIGE	2-0	76	3/8 CIRCLE		30	-	12

 Round Bodied

 Cutting

 Taper Cut

Product Description

Composition Purified connective tissue (mostly collagen) derived from sub-mucosal layer of ovine intestine. Chromic catgut is tanned by treating with chromium salts to prolong its resistance to absorption and to color the material to dark brown

Coating None

Sizes U.S.P. size 2 to 3-0

Color Yellow, Beige

Absorption Profile 90 to 100 days


Needle Availability Available as Pre-cut sutures and in wide variety of Round body, Taper cut and Cutting needles

Tensile Strength

Effective strength retention for Chromic Gut is 21 to 25 days and for Plain Gut it is 10 to 12 days

11% more knot pull tensile strength Vs. U.S.P. limit ensures excellent support throughout the critical wound healing period

Knot pull tensile strength (KgF)¹

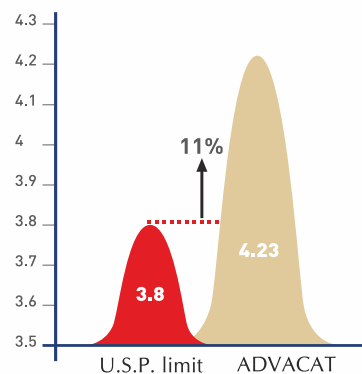

REFERENCES

¹: Data on file with AMS - Bench test done on chromic Gut U.S.P. size 1, 40mm, 1/2 Circle Round body



Advanced MedTech Solutions Pvt. Ltd.

Ceejay House, 305, Level 3, Shiv Sagar Estate, Dr. Annie Besant Road, Worli, Mumbai 400018. India.
T : (91)22-24946100 / 01 | F : (91)22-24946102 | W: www.amsltd.com | E : info@amsltd.com